


## FastGene® Culture Chamber (CC-01)

Protected environment always and everywhere

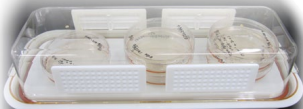

- ✓ **Especially designed to create a filtered environment for Petri dishes, 6-, 12-, 24- and 96-well plates**
- ✓ **The growth conditions are not altered**  
*The FastGene® Culture Chamber is equipped with exchangeable filters\* (CC-01F, pore size <math><0,4\mu\text{m}</math>), which are gas and humidity permeable. Hence, temperature, humidity, and  $\text{CO}_2$ -concentration are the same as the environmental. It even buffers in case of a quick opening of the incubator.  
*\*Filters should be changed every 2-3 weeks or after an infection in the incubator.**
- ✓ **Autoclavable for cell culture grade sterility**  
*The FastGene® Culture Chamber is autoclavable (lid, rubber ring and plate). Please sterilize when changing the filter sets.*